


# MicroSnap® | Coliform

MicroSnap® Coliform rapidly detects and enumerates Coliform bacteria. Single figure organisms can be detected in 8 hours, delivering results in the same working day or shift.


## Benefits:

- Get results in as few as 6 hours
- Large Enumerable Range
- Real-time data sync and alerts for early risk mitigation
- Data integration with EnSURE Touch and SureTrend
- Easy sample collection on liquid, solid and surface tests
- Snap Valve Technology
- Compact product design taking up minimal lab space
- Reduced cost compared to Petrifilm and third-party labs
- Biochemical detection
- Read using the highly sensitive EnSURE Touch luminometer
- Optimized workflow with the fewest hands-on steps
- AOAC Research Institute *Performance Tested Method*<sup>SM</sup> (#071302) validated for raw ground beef, RTE ham\*, cooked and raw chicken, raw cod, raw prawns, bagged lettuce\*, raw milk, sandwich, bottled water.
- \*Certified for enumeration only
- Store refrigerated (2 to 8 °C)

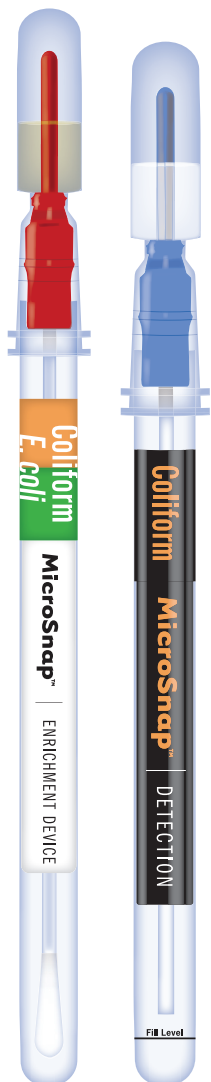


## Sensitivity & Specificity:

- Detects down to <10 CFU in 8 hours
- Detects Coliform bacteria including *Escherichia*, *Klebsiella*, *Citrobacter*, and *Enterobacter*

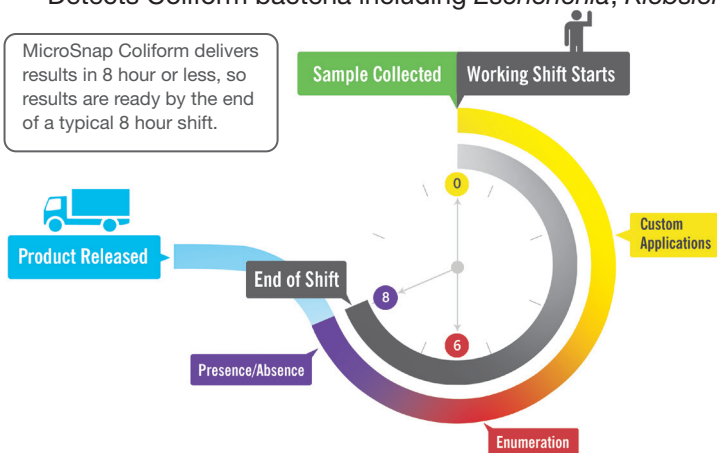

Catalog No.	Description	Quantity
MS1-CEC	MicroSnap® Coliform & <i>E. coli</i> Enrichment Device	100
MS1-EB-BROTH-9ML	MicroSnap® Enhanced EB Broth	25
MS2-COLIFORM	MicroSnap® Coliform Detection Device	100

## MicroSnap Coliform Procedure

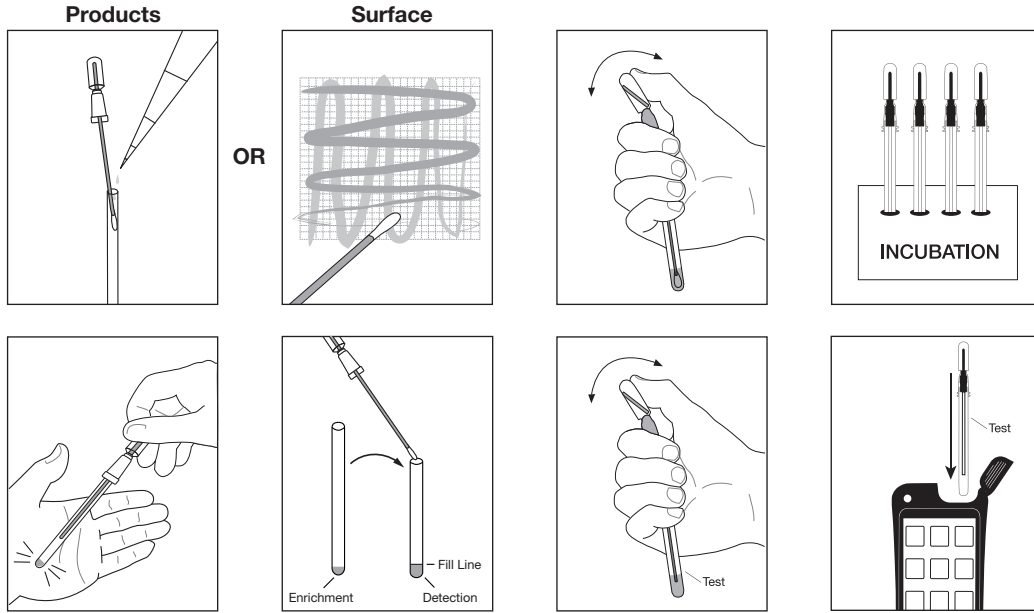

See an instructional demo at [HygienaTV](#) or [HygienaTV](#)

### Related Products

#### EnSURE® Touch Monitoring System

EnSURE® Touch is a high-sensitivity luminometer that measures multiple quality tests such as ATP, EB, Total Viable Counts, Coliform, *E. coli*, and more. Plus, send test results data to SureTrend® software to record, trend and analyze results over time. Get real-time notification and take your environmental monitoring program to the next level.



EnSURE® Touch Charging Dock

#### MicroSnap® *E. coli*

MicroSnap® *E. coli* detects and enumerates *E. coli* bacteria. Results are available in 8 hours or less. MicroSnap *E. coli* detects down to 1-10 CFUs.



#### MicroSnap® Total

MicroSnap® Total delivers total viable organisms test results in 7 hours. MicroSnap Total detects down to <10 CFU of gram-positive and gram-negative aerobic and facultative anaerobic bacteria.



#### MicroSnap® EB (*Enterobacteriaceae*)

MicroSnap® EB delivers test results in as few as 6 hours. MicroSnap EB detects down to 50 CFUs.



Catalog No.	Description	Quantity
ETOUCH	EnSURE® Touch Monitoring System	1
<a href="http://www.hygiene.com/suretrend">www.hygiene.com/suretrend</a>	SureTrend®	Inquire
MS1-CEC	MicroSnap® Coliform & <i>E. coli</i> Enrichment Device	100
MS2-COLIFORM	MicroSnap® Coliform	100
MS2-ECOLI	MicroSnap® <i>E. coli</i>	100
MS2-TOTAL	MicroSnap® Total	100
MS2-EB	MicroSnap® EB	100
CAL	CalCheck LED Calibration Verification Device	1